

Continuous Flow Microliter Control



- Excellent flow rate accuracy and reproducibility
- Continuous, constant flow with no pulsation for stable liquid transfer
- Flow rate adjustment between 5 and 5000 μL/min
- Maximum discharge pressure of 3 MPa (435 psi)
- No liquid leakage or exposure to outside air
- Transfer strong acids, alkalis, and water-based materials
- No change in flow rate even after long-term use
- Low dead volume (amount)
- Quick and easy maintenance for fast change-overs





- , .

_



- UP

Ε

S

[Control type]

Blank: Standard type
I: I/O signal control type

Model (flow rate standard)

5:5 mL/min

[Liquid-end part type]

KR : PEEK 6R : SUS316L [Power plug]

UP : Universal plug

[Language]

E: English Blank: Japanese [General]

S: Standard X: Special order

Pump specifications

Item		Q / QI-5	
Maximum flow rate*1		5000 μL/min (0.07 US G/h)	
Flow rate control range		5 to 5000 μL/min	
Maximum discharge pressure		3 MPa (435 psi)	
Connection type		1/4-28 UNF	
Connection diameter		φ1/16″×φ1/8″ ^{*2}	
Ambient temperature*3		0 to 40 °C	
Transfer liquid temperature		0 to 40 °C (no freezing)	
Transfer liquid viscosity		50 mPa⋅s or less	
Environmental resistance		IEC standard: Corresponds to IP65 (dust proofing and waterproofing)	
Power supply	Rated voltage	100 to 240 VAC±10%	
	Phase/frequency	1φ / 50 Hz or 60 Hz	
	Cord	1.8 m	
Weight		KR: 0.74 kg, 6R: 0.92 kg	
		·	

^{*1} Conditions : Clean water at room temperature. *2 Contact TACMINA for applications with different sizes of fittings and tubes.
*3 -10 to 50°C during transportation and storage.

Control functions

Function	Description		QI
Manual operation	Setting possible in units of 1 µL/min		
Auto-stop operation	Operating time: 1 to 9999 seconds or 1 to 9999 minutes Discharge volume: 1 to 9999 µL or 1 to 9999 mL		-
Pulse-input proportional operation	1 to 9999 µL/pulse or 1 to 9999 mL/pulse	_	•
Analog-input proportional operation	Controllable according to set target value (SV) and max. flow rate (HV)	_	

Liquid-end material

Item	KR	6R
Pump head	PEEK	SUS316L
Diaphragm	PTFE / CR	PTFE / CR
Check ball	Ruby	Ruby
Valve seat	Sapphire	Sapphire
Packing	PTFE	PTFE
Ball guide	PEEK	SUS316L

I/O signal specifications (QI only)


Signal		QI	
Analog	Input	1 port; 4 to 20 mA DC; Input resistance: Approx. 110 Ω	
Digital	Input	2 ports; No-voltage contact or open collector ^{*1} Maximum pulse count: 6000 pulses/min Minimum pulse width: 5 msec (ON time) Assigned to one of the following: Pulse, pause/run, level, or MAX operation	
	Output	2 ports; 25 VDC, 10 mA or less Assigned to one of the following: Unit pulse, alarm, or operation signal	


^{*1} Make sure the residual voltage of the open collector does not exceed 1 V.

Accessories

Tube (PTFE, φ1/16" × φ1/8") ····· 2 m	AC adapter (with A/C/O/BF plugs) ······· ×1
Syringe set ··································	Signal cable (QI only) · · · · 1.5 m
Joint assembly (ETFE) ······×2	Allen wrench (2.5 mm across flats) · · · · ×1
Nonslip sheet ···································	Operation manual ······ ×1

External dimensions

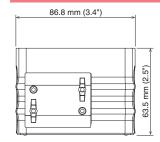

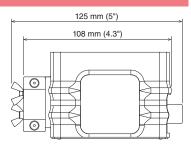

Consumables (per pump)

Part name	Quantity	Recommended replacement interval	
Diaphragm	2		
Check ball / valve seat assembly	4	4000 hours or 1 year	
Packing	8		

Wide variety of Milliliter Class pumps are also available. PVC type PTFE type Stainless steel type Milliliter Class pumps are able to handle flow rates between 0.1 and 100 mL/min. Choose from a wide range of liquid-end materials and control types to suit the application. Smoothflow Pump Q Series Milliliter Class

Product designs and specifications are subject to change without notice for product improvement.

TACMINA USA CORPORATION

105 W. Central Road, Schaumburg, IL 60195 U.S.A. Tel.844-822-6462 Fax.847-713-5410 URL www.tacminausa.com E-mail usasales@tacmina.com

TACMINA CORPORATION

Head Office:

2-2-14 Awajimachi, Chuo-ku, Osaka 541-0047 Japan Tel.+81-6-6208-3974 Fax.+81-6-6208-3978

URL www.tacmina.com E-mail trade@tacmina.com



